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OcClean Sonic Electric Toothbrush User Manual

Thank you for using OcClean Air2 Sonic Electric Toothbrush.

• Important Safety Tips

- 1 Due to the high-speed vibrations, it is normal to experience a tingling sensation while using the toothbrush for the first time. This feeling will gradually go away after 7 to 14 days.
- 2 If the toothbrush is not used for a long time, the battery will run out due to the power consumption from standby mode. This can be solved by charging, but the charging indicator on the brush will not light up instantly. Please wait for 6 - 10 minutes for the toothbrush to go back to the normal state.
- 3 It is normal for the toothbrush to heat up when it is charging. If the toothbrush becomes extremely hot, immediately stop charging and/or turn it off and send it to OcClean's authorized customer service center for checking.
- 4 The toothbrush can only be used with the OcClean special heads. OcClean will not be able to take any consequences caused by non-official brush head.
- 5 The toothbrush is for oral cleaning only. Please do not use it for other purposes.
- 6 The toothbrush is a personal cleaning item, please avoid sharing it with others.
- 7 Please store the toothbrush in a cool and dry place rather than anywhere exposes to direct sunlight.

Caution for Oral Sensitive Person:

- 1 The following 3 types of people are not recommended to use this toothbrush unless you have already consulted the professional dentists.
 - 2 The person who has undergone or is under treatment for periodontal surgery.
 - 3 The person with poor fitting dentures, continuously bleeding gums or any other oral related symptoms.
- The person whose oral is sensitive.

Unless directed by a guardian, otherwise this toothbrush is not suitable for people with reduced mental and physical capabilities, and those who lack related experience, and the children under the age of 12 are included.

Notes:

- * Please shut off the power when the charging base is needless or before cleaning.
- * Please regularly clean the power plug to prevent the accumulation of dust and dirt which can reduce the charging effect.

• Interaction Introduction

- 1 **Power On**
Start the toothbrush by short pressing the power button, the toothbrush will automatically shut down after a 2-mins brushing plan is finished.
- 2 **Intensity Switch**
OcClean Air2 has preset 2 types of 2mins brushing plan: Standard Cleaning(Default) and Deep Cleaning.
To switch the intensity, just long press the power button for 2s under working condition.
- 3 **Battery Indication**
Please charge the toothbrush when the red indicator light is on.

Bristle

Brush Head

ON/OFF Button

Short Press: Power On/Off
Long Press: Switch brushing intensity within standard mode and deep mode

Toothbrush Handle

There is a gap of about 0.5mm between the brush head and the handle to reserve a reasonable working space for the brush head.

Power Indicator

Charging status :
The white light flashes
Full Charge:
The white light shines
Low Battery:
The red light flashes fast


Annular Charging Interface

Magnetic Charging enables the toothbrush to be attracted and charged when it is close to the charging base.
Note: Please turn off the power when the charging is completed.


• Recommended Tooth Brushing Method

Brushing Time

Oral experts (Dentists) recommended that each brushing should be 2 to 3 minutes long.

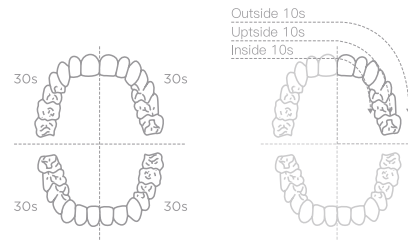• Recommended Gestures for Different Tooth Region

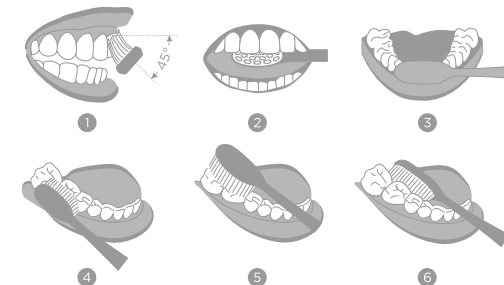• Using Instructions

Firstly, please insert the brush head forwardly onto the toothbrush body well. A small gap of about 0.5mm should be left between the head and the main body so the brush head can function properly.

Secondly, please turn on the toothbrush when it has been put into your mouth with a small amount of soaked toothpaste to prevent toothpaste from splashing.

Thirdly, please make sure the bristles at 45° from the gum line and move the toothbrush slowly.

Fourthly, hold the toothbrush handle and slowly move it in a circular motion. Do not press too hard.

Finally, please gargle with clean water after brushing.

● Specifications	
Product Name	Oclean Air2 Sonic Electric Toothbrush
Model Name	Oclean Air2
Input Power	2.5W
Charging	Oclean USB Charger
Input	5V ≈ 0.5A
Charging Standard	DC 5V 1A or above
Battery Life	30 Days (*)
Battery Capacity	800 mAh (**))
Charging Time	2.5 Hours (***)
Standard	GB 4706.1-2005 GB4706.59-2008

* Validated by 2 minutes brushing period, twice a day with low intensity mode.

** The battery capacity theoretically has a capacity range of 350mAh.

*** It is suggested to use a USB charger of 5V ≈ 0.5A.

● Charging Process

0%-95%, the power indicator light flashes.
95%-100%, the Power Indicator light is constantly on.

First Use

Charging time depends on the battery life demand of user.

Working Process


If battery power is less than 5%, the red indicator light will flash as a charging reminder.

Notes:

The battery cannot be replaced by users themselves and can only be replaced by the staff in Oclean Authorized Client Center.
Unauthorized bodies are forbidden to disassemble this product. The USB cable cannot be replaced by unauthorized body. If the product was damaged, please contact the Oclean After-sale Center for maintenance.


Please use the USB charger of 5V/1A or 5V/2A of mobile phone to charge. Such as the iPhone Charger, Tablet Charger etc.
If the charge current is less than 1A, it may prolong the charging time.
If the charge voltage is higher than 5V, it may damage the charging circuit.

● Cleaning and Maintenance

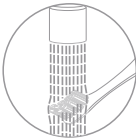
Cleaning The Brush Head

Clean the brush head after use and air dry it.
Do not pull the bristles, because it may result in deformation of the brush head.

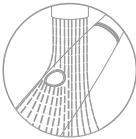Cleaning The Body

Use clean water to wash the surface of handle and wipe it with soft cloth and air dry it. Do not press the rubber sealing ring with sharp item or metal object, otherwise it will be damaged.

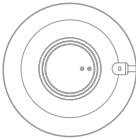Cleaning The Charger

Wipe the dirty parts of the charger with a soft cloth. Pull out the battery plug from the socket prior to cleaning the charger.
Wipe the surface of charger with wet cloth.

Changing The Head

It is suggested by dentists that the head should be changed every 2-3 months to maintain optimal brushing.

● Common Problems and Solutions

- The following are some common problems the users can face while using the toothbrush.
These problems can be easily solved by the users with the solutions given below. Please conduct the following trouble shooting techniques prior to contacting professional after-sale services.

The vibration strength becomes weak

Possible Problems:

- * The toothbrush is pressed too hard against the tooth.
- * The gap between the brush head and body is too loose or too tight.

Solutions:

- ✓ Lightly place the brush head on tooth. Pull the brush head from the body or gently insert the brush head leaving enough space between the brush head and body.

- The toothbrush functioning for short time

Possible Problems:

- * The charging cradle is not in contacted with the charging ring of the body, and may have resulted in improper charging which can also be detected by the LED light.
- * Insufficient battery charge.

Solutions:

- ✓ Clean the charger cradle.
 - ✓ Charge for longer period.
- Battery has completed its service life. If the red light of the foundation is on, please contact the Oclean Official Authorized Center for a replacement charger.

- The product is unable to work

Possible Problem:

- * Insufficient battery power.

Solutions:

- ✓ Please charge the toothbrush
- ✓ Please contact Oclean Official Support Center for a replacement battery if the first way doesn't help.

For detailed consultations,
Email: service@oclean.com

● FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

Guarantee Certificate

User Information

Name: _____ Tel: _____ E-Mail: _____
Address: _____ Zip Code: _____

Sales information

Name of Retailer: _____ Tel: _____ Sales Date: _____
Address: _____ Zip Code: _____ Invoice No: _____

Product Information

Type: _____ Production Place: _____

Warranty Card

Delivery Date	Fault description/ Reason/Disposal	Employed equipment	Acceptance date	Serviceman

Please ask for a maintenance report from the maintenance center and properly keep it with the warranty certificate. Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers guarantee for customers who purchase Oclean Smart Sonic Electric Toothbrush. In case of inconsistency between the service articles and the provisions of guarantee, the provisions of guarantee and the contents in the guarantee certificate shall govern.

Shenzhen Yunding Information Technology Co., Ltd.

32G, Tower 3, Dachong Business Building(The Third Phase), Shennan Avenue North, Tong Gu Road West, Yuehai Street Communities, Nanshan District, Shenzhen, Guangdong, China Zip Code: 518057
If there is fault occurring in the product you purchase, please contact Email: service@oclean.com

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Remark

○ The content of toxic and harmful substance in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011;

x The content of toxic and harmful substance in at least one of the homogeneous materials of this component is over the limited quantity as specified in GB/T 26572-2011;

Most materials of this product adopt environment-friendly materials. You can use this product without worry of infringing environment-protection policies.