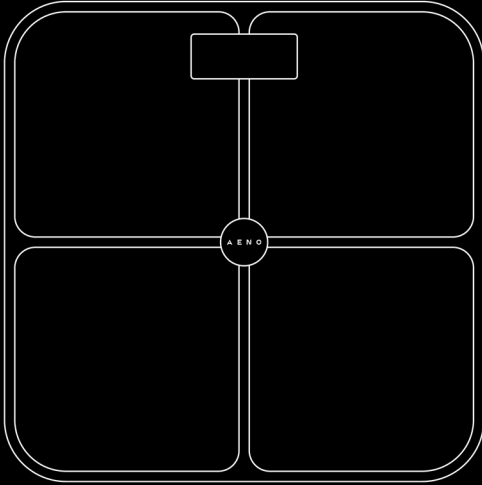
A E N O

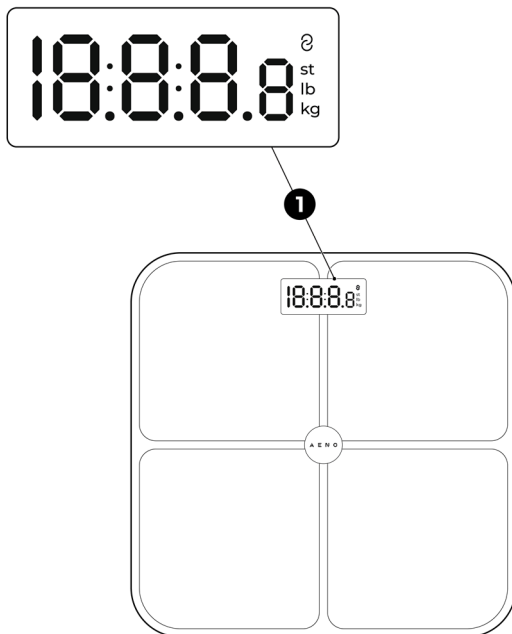
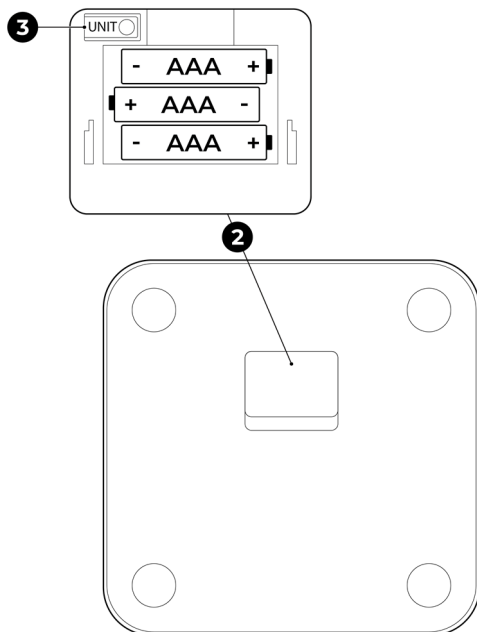
Quick Start Guide

Smart Body Scales BS1S

(Model: ABS0001S)

ENG	4	KAZ	18
ARA	5	LAV	19
BOS	6	LIT	20
BUL	7	NLD	21
CES	8	POL	22
DEU	9	POR	23
ELL	10	RON	24
EST	11	RUS	25
FRA	12	SLK	26
HRV	13	SLV	27
HUN	14	SPA	28
HYE	15	SRP	29
ITA	16	UKR	30
KAT	17	UZB	31

A**B**

ENGAENOTM ABS0001S Smart Body Scales are designed to determine the user's weight.

Technical Specifications

Power supply: AAA batteries (3 pcs), 4.5 V, 0.2 W (max). BT 5.0. Operating frequency: 2,400–2,480 MHz. Transmitter power: +4 dBm / 2.5 mW. Measuring range: 5–180 kg. Weight increment: 0.1 kg. Measurement error: ±0.3 kg (for 0–50 kg); ±0.43 kg (for 50–100 kg); ±0.5 kg (for more than 100 kg). Units of measurement: kg, lb, st:lb. Display: VA; 63×27 mm. Automatic activation/deactivation. Indication of the BT connection to the app. Control: manual, remote (AENO app). Color: White. Casing material: HIPS plastic, tempered glass with ITO coating. Dimensions (L×W×H): 280×280×27.5 mm. Weight: 1.6 kg. Operating conditions: temperature +5...+40 °C, RH 40–85 % (non-condensing). Storage conditions: temperature –20...+60 °C, RH up to 90 % (non-condensing).

Scope of Supply

Body Scales **ABS0001S**, AAA batteries (3 pcs), Quick Start Guide.

Restrictions and Warnings

Handle the glass surface of the scales with care: avoid damaging, cracking and chipping the surface as broken glass can cause injury. Supervise the use of the device by children and persons with physical and mental disabilities. Do not allow children to play with the device. Place the scales on a hard, flat, horizontal surface to avoid inaccurate measurements. Do not place the scales on carpets or any other soft surface. Do not load the device beyond the maximum value of the measuring range. If the device is not used for an extended period of time, remove the batteries and store them out of the reach of children and animals. Do not use chemical and harsh detergents, abrasive pastes, products containing acids and solvents, or metal sponges to clean the device.

ATTENTION! Do not expose the device to water or other liquids. Do not immerse the device in water. Do not wash the device or any part of it under running water or in a dishwasher.

Device Elements (see fig. A, B)

1 – display; 2 – battery compartment; 3 – "UNIT" button.

Device Indicators

 Active pairing with the AENO app	 Weight display scale	 Selected weight measurement unit
--	--	--

Preparing Device for Use

Unpack the device and clean it (see "Cleaning and Maintenance"). Open the battery compartment (2) cover and insert the batteries into the compartment, observing the correct polarity. Press the "UNIT" button (3) to select the unit of weight measurement. After each press, the display (1) will show the current unit of measurement: "kg" – kilogram, "lb" – pound, "st:lb" – stone. Close the cover of the compartment. Place the device on a hard, horizontal surface.

Note. The selected settings are saved for current use and for the next time the device is turned on. To reset the settings, press and hold the "UNIT" button for 5 seconds. The display will show the code "CLR". Also, the settings are automatically reset when the batteries are removed.

Operation of the Device

Stand on the surface of the device. The scales will turn on automatically, and the display will show your weight in the selected unit of measurement. The scales will turn off automatically within 15 seconds after the measurement is complete.

ATTENTION! Information shown by the scales cannot be considered medical information.

Cleaning and Maintenance

Note. Remove the batteries before cleaning the device.

Wipe the device with a soft cloth dampened with clean water and then wipe dry. If the scales are very dirty, wipe them with a soft cloth dampened with a mild soap solution. Then wipe with a soft cloth dampened with clean water and wipe dry.

Connecting to AENO Mobile App

Download the AENO mobile app from Google Play or App Store and register a new account (or log in to an existing one).

- Stand on the scales to turn them on. Place your smartphone as close to the device as possible.
- On the home screen of the app, tap "+" and then tap "Add Device".
- Note.** If you receive system prompts to turn on any features on your smartphone, follow the prompts.
- The device will be found automatically, and its image will appear in the "Autosearch" section. Tap the "Add" button.
- Wait until the message indicating that the device has been successfully added appears. Tap "Done".

Additional information about the user is available in the app: body fat percentage, muscle and bone mass, body water level, body mass index, etc. For more information on the features available in the app, refer to the complete user manual at aeno.com/documents.

Troubleshooting

- The device does not turn on. Possible cause: batteries are missing or discharged; polarity was not observed when inserting the batteries. Solution: insert charged batteries into the compartment, observing the polarity.
- The device does not connect to the app. Possible cause: BT is not enabled on the smartphone; the smartphone is too far away from the scales; scales settings failure. Solution: enable BT on your smartphone; place the smartphone no more than 5 m away from the scales; reset the scales to factory settings.

Error Codes Shown on Display

"Lo" – low power supply charge level; "Err" – maximum allowable weight exceeded; "C" – weighing failure.

ATTENTION! If none of the suggested methods have solved the problem, contact your supplier or an authorized service centre. Do not disassemble the device or attempt to repair it yourself.

If you have any questions or difficulties in using your AENO™ device, please e-mail support at support@aeno.com or chat with us at aeno.com/service-and-warranty. The specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

ASBISC reserves the right to modify the device and to make changes and additions to this document without prior notice to users.

The warranty period and service life is 2 years from the date of retail sale of the device.

Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Made in China. All trademarks and brands mentioned herein are the property of their respective owners. The date of manufacture is indicated on the package.

Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information about companies accepting quality claims and warranties is available at aeno.com/documents.

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