## CE Conformity Statements

This device complies with CE & R&TTE Directive. A copy of the EU Declaration of Conformity is available online at <http://www.picooc.com/en/suport.html>

*McGrp.Ru*

PICOOC Technology Co.,Ltd

PICOOC is a trademark of PICOOC Technology Ltd.

Bluetooth is a trademark of Bluetooth SIG,INC, Other product and Company names mentioned herein may be trademarks of their respective companies.

Made in PRC.designed by PICOOC c2014 All Rights Reserved.

Room 309, C building,No.28, Xingjiekou Outer Street,Xicheng District,Beijing,China

User Manual

Read the instruction manual carefully before using the device.

## CONTENTS

01-04	.....	English
05-08	.....	Deutsch
09-12	.....	Français
13-16	.....	Español
17-20	.....	Italiano
21-24	.....	Nederlands
25-28	.....	Polskie
29-32	.....	Русский
33-36	.....	Português
37-40	.....	Svenska

## Connect

## 01 Open the package and install the batteries

## Package List

PICOOC Smart Scale Mini


3 AAA 1.5V batteries

User manual

Remove the transparent protective film, open the battery cover under the scale, install the batteries and close the cover.
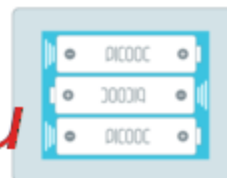
## 02 Scan the QR code to download PICOOC APP, register an account and then operate according to the guidance.

For old accounts, please select "Setting – Equipment Management" in PICOOC APP, chose the "+" symbol and then operate according to the guidance.

You may also search for PICOOC in APP Stores or Google Play.

\* Specific connecting methods may change due to APP upgrades, please refer to the latest version.

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