# BAMBA

## TOOTHCARE 1200 JET PRO

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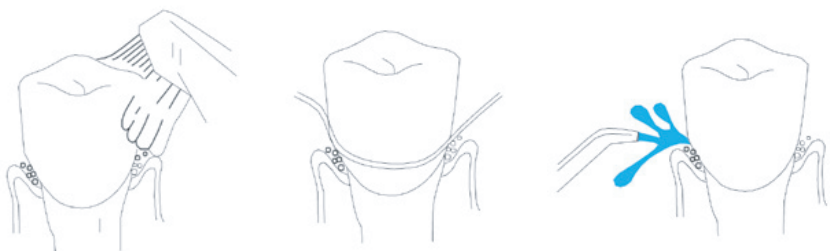
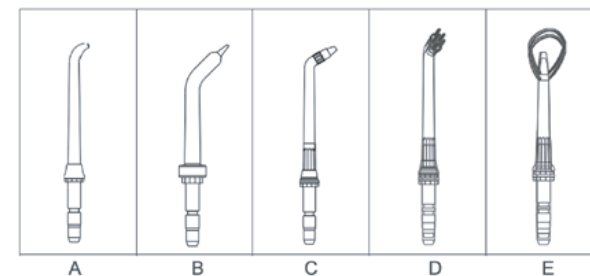
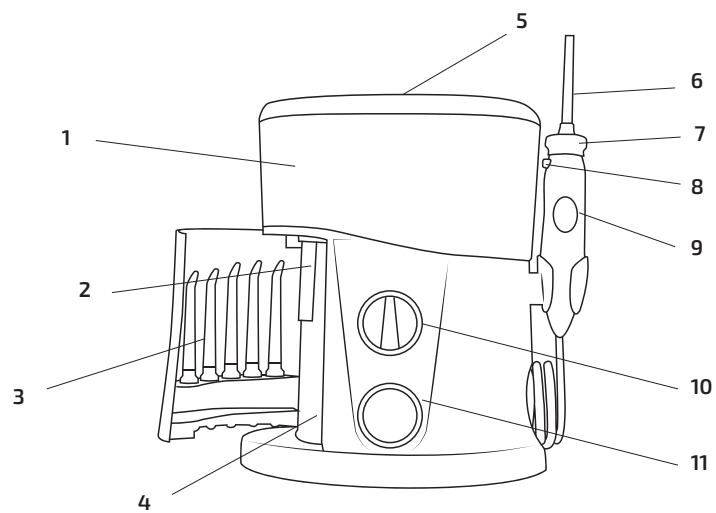
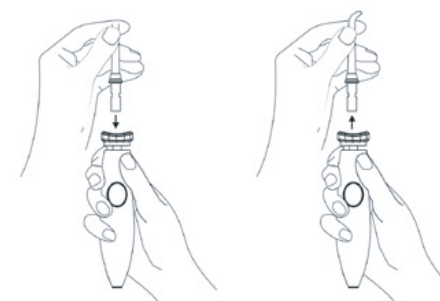
## SAFETY INSTRUCTIONS

**WARNING:** To avoid risk of fire, electric shock or injury:

- a. Do not look directly at the UV light.
- b. Do not let the device unattended while it is connected. Disconnect the plug when it is not being used.
- c. Do not use the appliance if its cord or plug are damaged, nor after it malfunctions or has been dropped or damaged in any way. Stop using the product immediately!

This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children must not play with the appliance. Cleaning and user maintenance must not be carried out by unsupervised children.

This device can only be used with cold or warm water or specific solutions.

Make sure that the mains voltage matches the voltage stated on the rating label of the appliance and that the wall outlet is grounded.

Check the power cord regularly for visible damage. If the cord is damaged, it must be replaced by the official Technical Support Service of Cecotec in order to avoid any type of danger.

Do not operate the product if its cord, plug or body show visible damage, do not operate properly or have been dropped.

Always turn off and unplug the appliance when not in use, before moving it to a new location and before assembling or disassembling it and before cleaning it. To avoid tipping over, wind the cord safely.

Do not direct water under the tongue, ear, nose or other sensitive areas. This product can produce pressures which may cause serious damage in these areas. Read this instruction manual for a correct operation.

Do not use it if you are wearing oral jewellery, remove it before using it.

Do not use it if you have an open wound in your mouth or tongue. If you have been advised by your physician to receive antibiotic pre-medication prior to dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.

Do not pull the cord nor use it for transport, do not use it as a handle, do not close a door over the cord nor pull it around edges or sharp corners. Do not place the device over the cord. Keep the cord away from heat sources.

Pull from the plug, not the cord to unplug it. Do not handle the cord or device with wet hands.

Turn off the product before unplugging it.

Unplug the irrigator before cleaning the device.

This device is indented for household use only. Do not operate outdoors. Do not use the irrigator in places where oxygen or aerosols are being used.

A normal phenomenon that causes bleeding gums is the use of a dental irrigator for the first time. It will depend on the gums different support capacity.

Store the appliance and its instruction manual in a safe and dry place when not in use.

Do not place or store the product in a place where it may fall or be dragged to a bathtub, shower or sink. Do not immerse in water or any other liquid. Do not use it while you shower.

Do not reach for an appliance that has fallen into water, unplug immediately.

Do not use accessories not provided by Cecotec. Use only the following accessories: high pressure normal nozzle, priodontal nozzle, braces nozzle, dental plaque nozzle and tongue cleaner nozzle. Fil the tank only with water or other recommended oral solutions.

### **Traditional cleaning**

Traditional toothbrushing cannot easily reach the unseen food remains and harmful bacteria. Fig. 1  
Floss has limited access below the gumline. Fig. 2

### **Oral Irrigator**

Remove trapped food debris and harmful bacteria with the irrigator. Fig. 3  
When using electric products, specially in the presence of children, you must follow basic safety precautions.

## 1. PARTS AND COMPONENTS

Fig. 4

1. Tank
2. Disinfection UV light
3. 7 nozzles
4. Safety lock
5. Tank cover
6. 360° rotary nozzle
7. Handle knob
8. Nozzle release button
9. Pause button
10. UV/ON/OFF knob
11. Pressure control knob with 10 pressure settings

Fig. 5

- A. High pressure regular nozzle x3
- B. Periodontal nozzle x1
- C. Braces nozzle x1
- D. Dental plaque nozzle x1
- E. Tongue cleaner nozzle x1

## 2. BEFORE USE

- Take the device out of its package and remove all packaging material.
- Check carefully if there is any damage. If there is, immediately notify the Technical Support Service of Cecotec.

### Box contents:

- Bamba ToothCare 1200 Jet Pro oral irrigator
- 7 nozzles
- Instruction manual

## 3. OPERATION

### UV disinfection

Use the disinfection UV light to clean the nozzles after every use. This system provides a controlled disinfection.

It has been checked that removes 99.9% of bacteria and viruses.

To activate the disinfection system, close the nozzle compartment door and turn the UV/ON/OFF knob to the UV position.

The UV light will only turn on if the door is well closed. The light will turn off if the door is open during the disinfection cycle. When closing it again, the cycle will start again.

The disinfection cycle acts for 15 minutes and then it turns off automatically.

Stop using the disinfection cycle if the UV light is on when the door is open or if the door is broken. The UV light can be harmful to the human eye and skin.

### Preparing the tank

Remove the tank from the base. Open the cover and fill the water tank with warm water.

Place the tank again in the base, press firmly and close the cover.

The tank cover can be completely removed to make its cleaning easy.

### Inserting and removing nozzles Fig. 6

Insert the nozzle in the centre of the knob in the upper part of the handle and press until it clicks. When it clicks you will see the knob lowers and stays fixed in this position. The color ring will be almost touching the end of the knob if the nozzle is correctly blocked in place.

To remove the nozzle from the handle, press the nozzle release button and remove the handle knob nozzle.

### Pressure setting

Turn the pressure control knob to the lowest position for the first use. Increase pressure gradually over time until the desired setting, or as directed by your dental professional.

### Use mode

Select the nozzle you want to use and insert it as previously indicated.

Connect the plug to a power supply.

Lean over the sink and place the nozzle on your mouth. Slightly close your lips to avoid splashes, but allow water to flow from your mouth to the sink.

Aim the nozzle at your teeth and with your free hand, turn the UV/OFF/ON switch to the ON position.

Set the pressure with the pressure control knob and the jet direction by turning the knob on the handle. Fig. 7

Direct the water jet in a 90° angle towards the gum.

To achieve better results, start in the molar area towards the front teeth. Fig. 8

Slide the nozzle over the gumline and stop between your teeth.

Keep the process until every area around and between teeth is clean.

**Pause/Control**

You can stop water temporarily at any time by pressing the Pause button on the handle. Fig. 9  
When finished, you must turn off the unit by placing the UV/OFF/ON knob in the OFF position and emptying the remaining liquid in the tank. Unplug the unit.

**Periodontal nozzle** Fig. 10

The periodontal nozzle releases water and anti-bacterial solution in the periodontal pockets. It is designed to reach the deepest areas in the periodontal cavity.

To use this nozzle, lean over the sink and place the nozzle carefully under the gumline, in a 45° angle. Turn on the unit with the lowest pressure setting and clean the gumline by introducing the nozzle between your teeth. Continue over the gumline until you finish cleaning.

**Braces nozzle** Fig. 11

The special nozzle for braces is exclusively designed to brush and rinse simultaneously the areas difficult to access around braces and other dental works.

To use the braces nozzle, lean over the sink and use the bristles to slightly clean around the braces, water will flow through the brush to rinse.

**Tongue cleaner nozzle** Fig. 12

Lean over the sink, place this nozzle over your tongue and turn on the unit with the lowest pressure setting. Increase it if necessary. Place the nozzle on the centre of your tongue and drag it forward applying slight pressure.

**Dental plaque nozzle** Fig. 13

Direct the jet in a 90° angle to the gumline, slightly close your lips to avoid splashes, but allow water to flow freely from your mouth to the sink. Start in the molar area (back teeth) and continue towards the front teeth, place the nozzle close to your teeth so bristles touch them softly.

Gently slide the nozzle along the gumline, passing briefly between the teeth to brush them gently and allow the water to flow between your teeth. Continue until all areas around and between the teeth have been cleaned, both inside and outside, top and bottom.

**Using other solutions**

You can use mouthwashes and antibacterial solutions intended for this purpose, paying special attention when cleaning.

To avoid clogging the irrigator after the use of solutions you should rinse it as follows: fill the tank partially with warm water and operate the device with the nozzle pointing towards the sink until the tank is empty.

**To finish**

Turn off the irrigator and completely empty the liquid left in the tank.  
If you use special solutions, clean and rinse thoroughly the tank.

**4. CLEANING AND MAINTENANCE****Cleaning**

Before cleaning it, turn off and unplug the irrigator.

After using special solutions in the tank, rinse the unit thoroughly.

After every use, you can use the UV light to sanitise the nozzles as explained above.

Fill half the tank with warm tap water, with the nozzle pointing towards the sink, turn the unit on and let it run until the tank is empty, this will help prevent the nozzle from clogging and help your unit last longer.

Clean the product, if necessary, with a soft cloth and a soft non-abrasive cleaning agent.

Disconnect the appliance from the power supply before cleaning.

**Maintenance and repair**

After cleaning the oral irrigator, place the tank upside down on a flat surface to make sure it completely dries. When it is completely dry, place it again on the body, also clean and dry.

You must place it in the original box or a similar one when it is not going to be used for a long time.

If the robot does not work properly, contact the Cecotec Technical Support Service.

Do not try to disassemble or repair the product. If this causes any problems with the product, the warranty will not cover it.

**5. TROUBLESHOOTING**

Problem	Cause	Solution
Leak between the nozzle and the handle knob	The nozzle is not correctly assembled	Remove the nozzle and try placing it back again until it clicks
Periodontal nozzle tear	The nozzle is worn	Replace it with a new nozzle
The oral irrigator base hose handle leaks	Damaged hose	Contact the official Technical Support Service of Cecotec to replace it
The tank leaks	The water tank valve is missing or upside down	Replace the tank valve Turn the valve
Inadequate pressure	The tank is not assembled correctly. The nozzle is clogged	Make sure the tank is full of water and correctly placed Clean the nozzle following the instructions
The unit does not release water	The motor is not working	Make sure the power supply works, check it with another device (check the circuit breaker and/or wall switch)

## 6. TECHNICAL SPECIFICATIONS

Model: Bamba ToothCare 1200 Jet Pro  
 Product reference: 04315  
 Power: 18 W  
 Voltage and frequency: 100-240 V~, 50/60 Hz  
 IPX4  
 Made in China | Designed in Spain

## 7. DISPOSAL OF OLD ELECTRICAL APPLIANCES

The European directive 2012/19/EU on Waste Electrical and Electronic Equipment (WEEE) specifies that old household electrical appliances must not be disposed of with the normal unsorted municipal waste.

Old appliances must be collected separately, in order to optimise the recovery and recycling of the materials they contain and reduce the impact on human health and the environment. The crossed out "wheeled bin" symbol on the

product reminds you of your obligation to dispose of the appliance correctly.

If the product has a built-in battery or uses batteries, they should be removed from the appliance and disposed of appropriately. Consumers must contact their local authorities or retailer for information concerning the correct disposal of old appliances and/or their batteries.

## 8. TECHNICAL SUPPORT AND WARRANTY

This product is under warranty for 2 years from the date of purchase, as long as the proof of purchase is submitted, the product is in perfect physical condition, and it has been given proper use, as explained in this instruction manual.

The warranty will not cover the following situations:


- The product has been used for purposes other than those intended for it, misused, beaten, exposed to moisture, immersed in liquid or corrosive substances, as well as any other fault attributable to the customer.
- If the product has been disassembled, modified, or repaired by persons not authorised by the Cecotec official Technical Support Service.
- Faults deriving from the normal wear and tear of its parts, due to use.

The warranty service covers every manufacturing defect of your appliance for 2 years, based on current legislation, except consumable parts. In the event of misuse, the warranty will not apply.


If at any moment you detect any problem with your product or have any doubt, do not hesitate to contact Cecotec official Technical Support Service at +34 963 210 728.

## 9. COPYRIGHT


The intellectual property rights over the texts in this manual belong to CECOTEC INNOVACIONES, S.L. All rights reserved. The contents of this publication may not, in whole or in part, be reproduced, stored in a retrieval system, transmitted or distributed by any means (electronic, mechanical, photocopying, recording or similar) without the prior authorization of CECOTEC INNOVACIONES, S.L.

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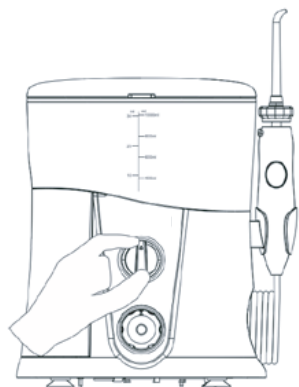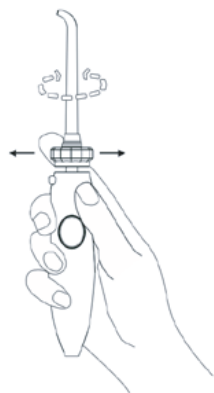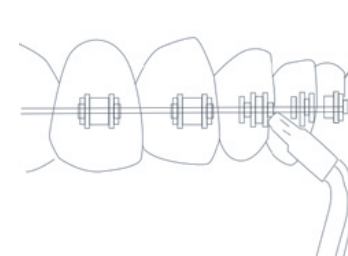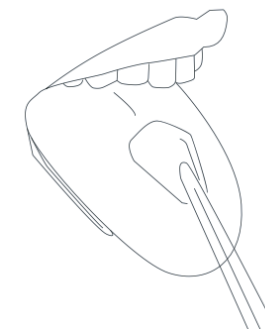
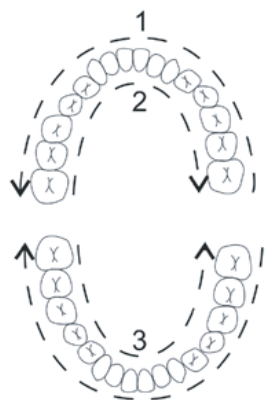
[www.cecotec.es](http://www.cecotec.es)

Cecotec Innovaciones S.L.  
C/ de la Pinadeta s/n, 46930  
Quart de Poblet, Valencia (Spain)  
PE01210924