
YUNMAI Massage Gun mini 2

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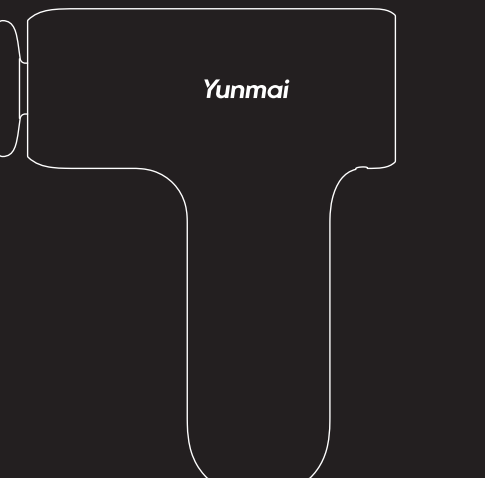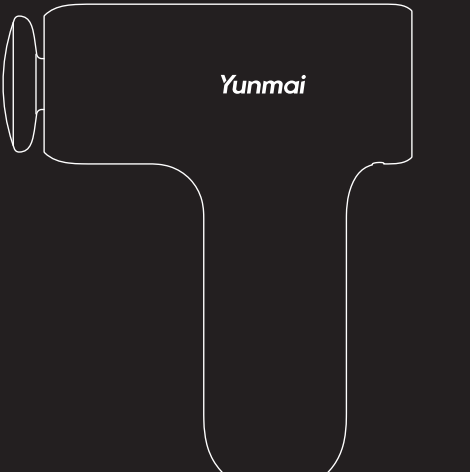
Instagram: [yunmaiglobal](https://www.instagram.com/yunmaiglobal)

Facebook: Yunmai

Manufacturer: Shenzhen Yunmai Technology Co., Ltd.

Address: Room 01, 18-19 / F, Chiwan Headquarters Building
Zhaoshang Street, Nanshan District, Shenzhen 518000
Guangdong P.R. China

Please read this manual carefully before use. For future reference, please ensure this manual is properly kept.

This device is designed for intensive percussion (shock) massage of various parts of the human body, helps overcome fatigue, tension, stress, can increase blood circulation, speed up the recovery process after training or a hard day at work, increase range of motion and muscle elasticity, eliminate negative sensations and discomfort.

Key Information

To reduce risks of malfunction, electric shock, injury, fire, death and damage to equipment, please read all statements and warnings in this manual carefully before use.

Prohibitions

- The use of this product by underaged and injured persons is strictly prohibited. In the event of muscle injury, please consult your physician before using this product.
- Do not use this product on wet or dirty skin. This product can be used on clothed areas. Users can press this product on the skin gently, move it along different spots, letting it stay at each particular spot for about 60 seconds.
- Do not use this product on body parts above the neck, nor on body parts with less muscle or bone content. This product can only be used on soft tissues when the user does not experience pain or adverse reactions.
- Do not use this product on body parts with bruises, pain or adverse reaction.
- Do not place other objects in the vent of the product.
- Do not immerse the product in water. Do not allow water to enter the vents of the product.
- Do not place the product in or near sources of fire.
- Do not throw or drop the product on purpose.
- Do not disassemble the product.
- Please do not put your fingers, hair or other parts of the body close to the gap between the massage head and the massage gun to avoid being trapped while the device is in operation.
- This product is intended for adults aged 18 and above. Minors and individuals over 60 years old should use it under supervision to prevent accidents. Children under the age of 8 are strictly prohibited from using this product.
- This product is not intended for independent use by individuals with limited physical or mental capabilities, or with insufficient experience and knowledge, unless they have been trained and instructed on the use of the product by individuals responsible for their safety.


Statement

Those having any one of the following circumstances shall consult the doctor before using the product.


- Pregnant women
- Pacemaker wearer
- Undergone surgery recently
- Abnormalities in spines, joints or bones
- Complications caused by other diseases (such as neuropathy or retinal damage)

Current research focuses on the influence of vibrational training on special medical conditions and may reduce the scope for the above-mentioned contraindications. Practical experience has shown that the combination of treatment and vibrational training also serves as a type of therapy in many circumstances. However, such an approach should be recommended by physicians, experts or therapists or be conducted in their presence.

Product Introduction

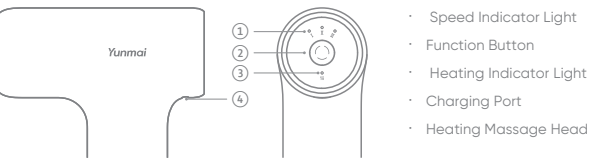
Packing List

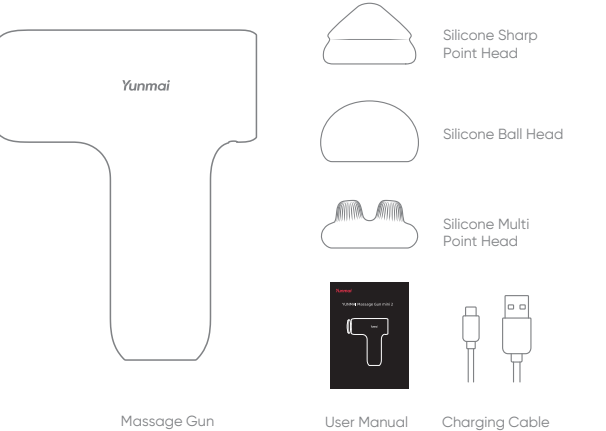Operation Guide


(1) Installation

Please install the silicone massage head when the product is turned off, as shown in the picture, rotate the silicone massage head onto the massage gun.

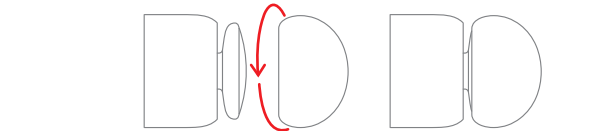
Warning

The heating massage head is not removable. Please install or remove the silicone massage head when the device is turned off. Ensure that the silicone massage head is installed firmly before starting the device. Do not adjust or remove the silicone massage head during use to avoid pinching your fingers.

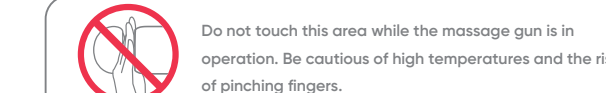This product has a heating surface, and it must be used with caution by individuals insensitive to heat.


(2) Power On, Speed Switch, Heating Switch, Power Off

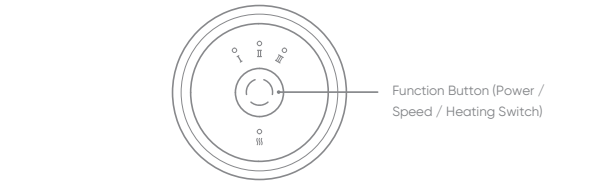Indicator Lights


Power On: Press and hold the function button for 2 seconds to turn on the device.
Speed Switch: Press the function button to switch the vibration speed. There are 3 speeds to choose from.
Heating Switch: Double-click the function button to switch to heating mode.
Power Off: After powering on, press and hold the function button for 2 seconds to turn it off. The device will automatically turn off after 2 minutes of no operation.

(3) Indicator Lights

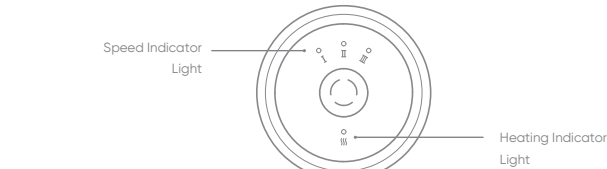Speed Indicator Light			Heating Indicator Light		
Color	Status	Description	Color	Status	Description
White	○ ● ●	Speed 1	White	○	Normal Mode
White	○ ○ ●	Speed 2	Red	●	Heating Mode
White	○ ○ ○	Speed 3			

(4) Charging

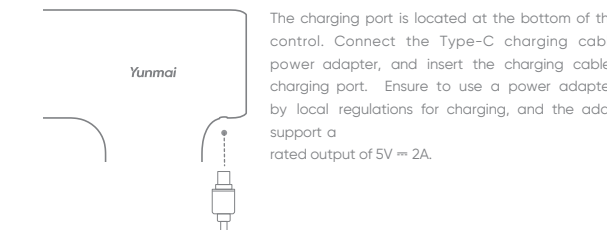Non-charging State			Charging Status		
Indicator Light color	Status	Battery Level	Indicator Light Status	Status	Battery Level
Red	● ● ●	<10%	White Light Blinking	○ ● ●	5%~49%
			White Light Blinking	○ ○ ●	50%~79%
			White Light Blinking	○ ○ ○	80%~99%
			White Light On	○ ○ ○	100%

Notice

It is prohibited to use the massage gun while charging! When the indicator light is stably on, charging is complete. Please unplug the charging cable promptly.

(5) Instructions

- Select and install the appropriate massage head. Turn on the massage gun and adjust the speed according to personal preference, it is recommended to start with a low speed. Gently apply the massage head to the muscle area that needs relaxation, moving it back and forth slowly. If there's a specific sore spot, apply appropriate pressure to the sore spot and hold it there for 15 to 30 seconds.
- Be mindful of the inherent impact force of the massage gun; avoid excessive pressure to prevent potential muscle and bone damage.
- It is recommended to use the massage gun on well-developed muscle areas and avoid staying in one spot for more than 60 seconds.
- The recommended massage duration is 10 minutes. If soreness persists, you can extend the massage time as needed.
- Keep individual sessions under 15 minutes and continuous use under 30 minutes. Strictly follow the instructions in this manual for proper product operation.

Specification

Product Name: YUNMAI Massage Gun mini 2
Model No.: YMFG-M356
Battery Type: Rechargeable lithium battery
Battery Capacity: 1500mAh
Rated Input: 5V = 2A
Charge Interface: Type-C
Working Temperature: 0~40°C
Working Humidity: 30~75%
Dimensions: 120.5 × 138.4 × 47 mm Weight: Approximately 383g

Frequently Asked Questions

Troubleshooting of Massage Gun Operation	
Fault description	Troubleshooting
The product cannot be started	Check whether the power is sufficient (if the product has not been used for a long time, please charge it for more than 4 hours)
The product cannot be charged	Check whether the charging connector and power supply are plugged in place
The product stops suddenly after 15 minutes of continuous operation	This phenomenon is normal as it is caused by the setting for long-term operation protection, which switches from operation to standby mode. And the product can be restarted through short pressing
The product stops working after the massage head encounters the phenomenon of locked-rotor	This phenomenon is normal as it is caused by the setting for abnormal operation protection, which switches from operation to standby mode. And the product can be restarted through short pressing.

Reasonable use of the product can effectively extend its service life.

Suggestions:

- Do not throw the product during use.
- Keep the maximum continuous service time lower than 30 minutes.
- Please charge the product every 2 months if it is not used for a long time.

Storage and Transportation

- Please store the massage gun in a cool, dry place away from direct sunlight and moisture.
- Please avoid storing the massage gun in environments outside 30%~90% humidity, temperatures above 40°C, or below 0°C. Keep it away from moisture, water sources, heat sources, and corrosive environments to prevent battery performance issues.
- When storing or carrying the massage gun outdoors, please place it in a portable storage bag.
- Please avoid impacting, pressuring, dropping, throwing, or scratching the massage gun. This product is suitable for transportation by car, train, ship, airplane, etc.
- During the transportation, please ensure the massage gun is securely positioned upright, avoiding compression, and preventing severe vibrations, collisions, rolling, and falling. If the massage gun has not been used for a long time, please recharge it for at least 4 hours before reuse.

WEEE Information

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points.

Battery Directive - Information for Users

In the European Union, this label indicates that the batteries in this product should be collected separately and not be disposed of with household waste. Substances in batteries can have a potential negative impact on health and environment and you have a role in recycling waste batteries thus contributing to the protection, preservation, and improvement of the quality of the environment. You should contact your local authority or retailer for details of the collection and recycling schemes.

EU Declaration of Conformity

We, Shenzhen Yunmai Technology Co., Ltd., hereby declare that this equipment is in compliance with the applicable Directives and European Norms, and amendments. The full text of the EU declaration of conformity is available at the following internet address: <http://www.yunmaiglobal.com>